



A **GUI-based** framework for **segmentation**, **tracking** and **cell cycle annotations** of microscopy imaging data. It includes three of the most accurate deep learning methods, [Cellpose](#), [YeaZ](#), and [YeastMate](#).

**Cell-ACDC** can load **2D**, **3D** (either single z-stacks or 2D images over time) and **4D** (3D z-stacks over time) images.

*Written in Python 3 by Francesco Padovani and Benedikt Mairhoermann.*

\*Tested on Windows 10 64 bit, macOS , and Linux Mint 20.1

## Installation

1. Download the latest release from [here](#).
2. If you don't already have Python or Anaconda, download, and install Miniconda for Python 3.8 [here](#). We recommend using Anaconda even if you already have Python.
3. Follow the instructions below specific to your OS

### Installing on Windows using conda

1. Unzip the latest release you downloaded before. For this example, I will assume it was unzipped into `C:\Users\Frank`
2. Open the Anaconda Prompt (you should be able to find it from the search bar)
3. Navigate to the folder where you unzipped Cell\_ACDC, (in this example it is `C:\Users\Frank\Cell-ACDC`) by typing `cd "C:\Users\Frank\Cell_ACDC"` Press "Enter" to confirm. Note that if you unzipped into a drive different from `C:\` you first need to change the drive letter in your terminal. To do so, type the letter first (e.g., `G:`) and then you can navigate with the `cd` command.

```
Anaconda Prompt (Miniconda3)
(base) C:\Windows>cd "C:\Users\Frank\Cell_ACDC"
(base) C:\Users\Frank\Cell_ACDC>
```

4. Now type the following commands **one at the time** (press "Enter" after each command and type "Y" when requested):

```
conda update conda
```

```
conda clean --all
```

```
conda env create --file acdc.yml
```


Anaconda will create the environment with Python 3.8 and all the packages required. This step can take several minutes (about 20 minutes if I have to guess, but it depends on your internet connection speed).

If successful, your terminal should now look like the screenshot below (red circle around the part that will tell you that the installation was successful). If you had an error, you could try installing using pip (see instructions below) or open an issue [here](#).

```
Anaconda Prompt (Miniconda3)
openssl-1.1.1l | 5.7 MB | ##### 100%
snappy-1.1.8 | 50 KB | ##### 100%
debugpy-1.4.1 | 3.2 MB | ##### 100%
pyzmq-22.2.1 | 464 KB | ##### 100%
h5py-3.4.0 | 1.1 MB | ##### 100%
winpty-0.4.3 | 1.1 MB | ##### 100%
Preparing transaction: done
Verifying transaction: done
Executing transaction: | "By downloading and using the CUDA Toolkit conda packages, you accept the terms and conditions of the CUDA End User License Agreement (EULA): https://docs.nvidia.com/cuda/eula/index.html"
done
Installing pip dependencies: \ Ran pip subprocess with arguments:
['C:\Users\Frank\conda\envs\acdc\python.exe', '-m', 'pip', 'install', '-U', '-r', 'C:\Users\Frank\Cell_ACDC\condaenv.0me8bvtc.requirements.txt']
Pip subprocess output:
Requirement already satisfied: PyQt5 in c:\users\frank\conda\envs\acdc\lib\site-packages (from -r C:\Users\Frank\Cell_ACDC\condaenv.0me8bvtc.requirements.txt (line 1)) (5.12.3)
Collecting PyQt5
  Using cached PyQt5-5.15.4-cp36.cp37.cp38.cp39-none-win_amd64.whl (6.8 MB)
Collecting opencv-python-headless
  Using cached opencv_python_headless-4.5.3.56-cp38-cp38-win_amd64.whl (34.8 MB)
Collecting PyQt5-sip<13,>=12.8
  Using cached PyQt5_sip-12.9.0-cp38-cp38-win_amd64.whl (63 kB)
Collecting PyQt5-Qt5>=5.15
  Using cached PyQt5-Qt5-5.15.2-py3-none-win_amd64.whl (50.1 MB)
Requirement already satisfied: numpy>=1.17.3 in c:\users\frank\conda\envs\acdc\lib\site-packages (from opencv-python-headless->-r C:\Users\Frank\Cell_ACDC\condaenv.0me8bvtc.requirements.txt (line 2)) (1.20.3)
Installing collected packages: PyQt5-sip, PyQt5-Qt5, PyQt5, opencv-python-headless
Attempting uninstall: PyQt5-sip
  Found existing installation: PyQt5-sip 4.19.18
  Uninstalling PyQt5-sip-4.19.18:
    Successfully uninstalled PyQt5-sip-4.19.18
Attempting uninstall: PyQt5
  Found existing installation: PyQt5 5.12.3
  Uninstalling PyQt5-5.12.3:
    Successfully uninstalled PyQt5-5.12.3
Successfully installed PyQt5-5.15.4 PyQt5-Qt5-5.15.2 PyQt5-sip-12.9.0 opencv-python-headless-4.5.3.56
done
#
# To activate this environment, use
#
# $ conda activate acdc
#
# To deactivate an active environment, use
#
# $ conda deactivate
#
(base) C:\Users\Frank\Cell_ACDC>
```

## Installing on Windows using pip

1. Download and install Python 3.8.4 from [here](#). **Make sure to check the option** Add Python 3.8 to PATH and then install with default options.

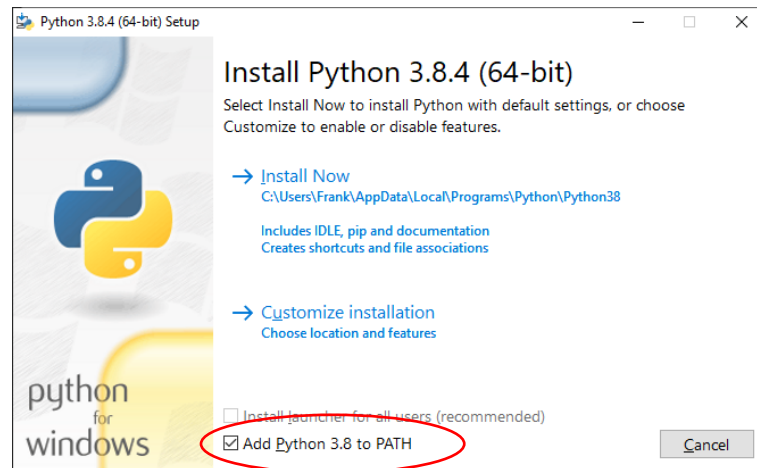

2. Unzip the latest release you downloaded before. For this example, I will assume it was unzipped into `C:\Users\Frank`
3. Open a terminal (either a Command Prompt or PowerShell, you can find both from the search bar)
4. Navigate to the folder where you unzipped Cell\_ACDC, (in this example it is `C:\Users\Frank\Cell_ACDC`) by typing `cd "C:\Users\Frank\Cell_ACDC"` Press *Enter* to confirm. Note that if you unzipped into a drive different from `C:`, you first need to change the drive letter in your terminal. To do so type the letter first (e.g., `G:`) and then you can navigate with the `cd` command.
5. Now type the following commands **one at the time** (press "Enter" after each command and type "Y" when requested):

```
py -m pip install --upgrade pip
```

```
py -m venv env
```


```
.\env\Scripts\activate
```

```
py -m pip install -r requirements.txt
```

You will now see all the required packages being installed. If successful, your terminal should look like the screenshot below (you can ignore that warning). If you had an error, you could try installing using Anaconda (see instructions above) or open an issue [here](#).

```
PowerShell
1 nbclassic-0.3.1 nbclient-0.5.4 nbconvert-6.1.0 nbformat-5.1.3 nest-asyncio-1.5.1 networkx-2.6.3 notebook-6.4.4 numba-0
.53.1 numexpr-2.7.3 numpy-1.19.5 oauthlib-3.1.1 opencv-python-4.5.2.54 opt-einsum-3.3.0 packaging-21.0 pandas-1.2.4 pand
ocfilters-1.5.0 parso-0.8.2 pickleshare-0.7.5 pillow-8.3.2 prometheus-client-0.11.0 prompt-toolkit-3.0.20 protobuf-3.18.
0 pyasn1-0.4.8 pyasn1-modules-0.2.8 pycparser-2.20 pyglet-1.5.16 pygments-2.10.0 pyparsing-2.4.7 pyqtgraph-0.12.1 pyrsis
tent-0.18.0 python-dateutil-2.8.2 pytz-2021.1 pywin32-301 pywinpty-1.1.4 pyzmq-22.2.1 requests-2.25.1 requests-oauthlib-
1.3.0 requests-unixsocket-0.2.0 rsa-4.7.2 scikit-image-0.18.1 scikit-learn-0.24.2 scipy-1.7.1 seaborn-0.11.1 six-1.15.0
sniffio-1.2.0 tables-3.6.1 tensorboard-2.6.0 tensorboard-data-server-0.6.1 tensorboard-plugin-wit-1.8.0 tensorflow-2.6.0
tensorflow-estimator-2.6.0 termcolor-1.1.0 terminado-0.12.1 testpath-0.5.0 threadpoolctl-2.2.0 tifffile-2021.8.30 torch
-1.9.0 tornado-6.1 tqdm-4.59.0 traitlets-5.1.0 typing-extensions-3.7.4.3 urllib3-1.26.6 wcwidth-0.2.5 webencodings-0.5.1
websocket-client-1.2.1 werkzeug-2.0.1 wheel-0.37.0 wrapt-1.12.1
WARNING: You are using pip version 20.1.1; however, version 21.2.4 is available.
You should consider upgrading via the 'C:\Users\Frank\Cell_ACDC\env\Scripts\python.exe -m pip install --upgrade pip' com
mand.
(env) PS C:\Users\Frank\Cell_ACDC>
```

## Installing on macOS using conda

1. Unzip the latest release you downloaded before. For this example, I will assume it was unzipped into `/Users/anikavanessaseel/Documents/GitHub/Cell_ACDC`
2. Open a **Terminal** (Click the Launchpad icon  in the Dock, type “Terminal” in the search field, then click Terminal)
3. Navigate to the folder where you unzipped Cell\_ACDC, (in this example it is `/Users/anikavanessaseel/Documents/GitHub/Cell_ACDC`) by typing `cd` plus drag-and-drop the `Cell_ACDC` folder from the Finder. Press *Enter* to confirm.
4. Now type the following commands **one at the time** (press “Enter” after each command and type “Y” when requested):

```
conda update conda
```


```
conda clean --all
```

```
conda env create --file acdc.yml
```

Anaconda will create the environment with Python 3.8 and all the packages required. This step can take several minutes (about 20 minutes if I have to guess, but it depends on your internet connection speed).

If you had an error, you could try installing using pip (see instructions below) or open an issue [here](#).

## Installing on macOS using pip

1. Download and install Python 3.8.5 from [here](#). Install with default options.
2. Unzip the latest release you downloaded before. For this example, I will assume it was unzipped into `/Users/anikavanessaseel/Documents/GitHub/Cell_ACDC`
3. Open a **Terminal** (Click the Launchpad icon  in the Dock, type “Terminal” in the search field, then click Terminal)
4. Navigate to the folder where you unzipped Cell\_ACDC, (in this example it is `/Users/anikavanessaseel/Documents/GitHub`) by typing `cd` plus drag-and-drop the `Cell_ACDC` folder from the Finder. Press *Enter* to confirm.
5. Now type the following commands **one at the time** (press “Enter” after each command and type “Y” when requested):

```
python3 -m pip install --user --upgrade pip
```

```
python3 -m venv env
```

```
source env/bin/activate
```

```
python3 -m pip install -r requirements.txt
```

You will now see all the required packages being installed. If successful, your terminal should look like the screenshot below. If you had an error, you could try installing using Anaconda (see instructions above) or open an issue [here](#).

```
-bash: src: command not found
Anikas-MacBook-Pro:~ anikavanessaseel$ cd documents
Anikas-MacBook-Pro:documents anikavanessaseel$ cd github
Anikas-MacBook-Pro:github anikavanessaseel$ cd yeast_acdc
Anikas-MacBook-Pro:yeast_acdc anikavanessaseel$ python3 -m venv env
Anikas-MacBook-Pro:yeast_acdc anikavanessaseel$ source env/bin/activate
(env) Anikas-MacBook-Pro:yeast_acdc anikavanessaseel$ python3 -m pip install -r requirements.txt
ERROR: Could not find a version that satisfies the requirement python==3.8.5 (from -r requirements.txt (line 1)) (from versions: none)
ERROR: No matching distribution found for python==3.8.5 (from -r requirements.txt (line 1))
WARNING: You are using pip version 20.1.1; however, version 21.1.3 is available.
You should consider upgrading via the '/Users/anikavanessaseel/Documents/GitHub/Yeast_ACDC/env/bin/python3 -m pip install --upgrade pip' command.
(env) Anikas-MacBook-Pro:yeast_acdc anikavanessaseel$ python3 -m pip install -r requirements.txt
Collecting pyqtgraph==0.12.1
  Downloading pyqtgraph-0.12.1-py3-none-any.whl (939 kB)
    ERROR: Could not find a version that satisfies the requirement pyqt==5.12.3 (from -r requirements.txt (line 2)) (from versions: none)
    ERROR: No matching distribution found for pyqt==5.12.3 (from -r requirements.txt (line 2))
    WARNING: You are using pip version 20.1.1; however, version 21.1.3 is available.
    You should consider upgrading via the '/Users/anikavanessaseel/Documents/GitHub/Yeast_ACDC/env/bin/python3 -m pip install --upgrade pip' command.
(env) Anikas-MacBook-Pro:yeast_acdc anikavanessaseel$ python3 -m pip install -r requirements.txt
Collecting pyqtgraph==0.12.1
  Using cached pyqtgraph-0.12.1-py3-none-any.whl (939 kB)
Collecting PyQt5==5.15.4
  Downloading PyQt5-5.15.4-cp38-cp38-macosx_10_13_intel.whl (7.0 MB)
    Collecting pandas==1.2.4
      Downloading pandas-1.2.4-cp38-cp38-macosx_10_9_x86_64.whl (10.5 MB)
    Collecting tables==3.6.1
      Downloading tables-3.6.1-cp38-cp38-macosx_10_9_x86_64.whl (4.3 MB)
    Collecting opencv-python==4.5.2.54
      Downloading opencv_python-4.5.2.54-cp38-cp38-macosx_10_15_x86_64.whl (43.7 MB)
    Collecting matplotlib==3.1.2
      Downloading matplotlib-3.1.2-cp38-cp38-macosx_10_9_x86_64.whl (13.2 MB)
    Collecting tqdm==4.59.0
      Downloading tqdm-4.59.0-py2.py3-none-any.whl (74 kB)
    Collecting natsort==7.1.1
      Downloading natsort-7.1.1-py3-none-any.whl (35 kB)
    Collecting scipy==1.6.2
      Downloading scipy-1.6.2-cp38-cp38-macosx_10_9_x86_64.whl (30.8 MB)
    Collecting scikit-image==0.18.1
      Downloading scikit_image-0.18.1-cp38-cp38-macosx_10_9_x86_64.whl (12.7 MB)
    Collecting pyglet==1.5.16
      Downloading pyglet-1.5.16-py3-none-any.whl (1.1 MB)
    Collecting seaborn==0.11.1
      Downloading seaborn-0.11.1-py3-none-any.whl (285 kB)
    Collecting numba==0.53.1
      Downloading numba-0.53.1-cp38-cp38-macosx_10_14_x86_64.whl (2.2 MB)
    Collecting requests==2.25.1
      Downloading requests-2.25.1-py2.py3-none-any.whl (61 kB)
    Collecting tensorflow==2.3.0
      Downloading tensorflow-2.5.0-cp38-cp38-macosx_10_11_x86_64.whl (195.7 MB)
```

## Installing on Linux using conda

1. Open a **Terminal**
2. Unzip the latest release you downloaded before. For this example, I will assume it was unzipped into `/home/elpado/GitHub/`
3. Navigate to the folder where you unzipped Cell\_ACDC, (in this example it is `/home/elpado/GitHub/`) by typing `cd "/home/elpado/GitHub/Cell_ACDC"`
4. Now type the following commands **one at the time** (press "Enter" after each command and type "Y" when requested):

```
conda update conda
```

```
conda clean --all
```

```
conda env create --file acdc.yml
```

Anaconda will create the environment with Python 3.8 and all the packages required. This step can take several minutes (about 20 minutes if I have to guess, but it depends on your internet connection speed).

If successful, your terminal should now look like the screenshot below (red circle around the part that will tell you that the installation was successful). If you had an error, you could try installing using pip (see instructions below) or open an issue [here](#).

## Installing on Linux using pip

1. Open a Terminal
2. Make sure you have Python 3.8 and pip installed. Check if you have Python with `python --version` command and check if you have pip with `pip help` command. If you don't have them install with the following commands:

```
sudo apt-get update
```

```
sudo apt-get install python3.8
```

```
sudo apt-get install python3-pip
```

3. Unzip the latest release you downloaded before. For this example, I will assume it was unzipped into `/home/elpado/GitHub/`
4. Navigate to the folder where you unzipped Cell\_ACDC, (in this example it is `/home/elpado/GitHub/`) by typing `cd "/home/elpado/GitHub/Cell_ACDC"`
5. Now type the following commands **one at the time** (press "Enter" after each command and type "Y" when requested):

```
python3 -m pip install --user --upgrade pip
```

```
python3 -m venv env
```

```
source env/bin/activate
```


```
python3 -m pip install -r requirements.txt
```

You will now see all the required packages being installed. If you had an error, you could try installing using Anaconda (see instructions above) or open an issue [here](#).

# First steps

## Starting the main launcher

1. Open a terminal:
  - Windows: **Anaconda Prompt** if you installed with conda
  - Windows: Command Prompt or **PowerShell** if you installed with pip
  - Unix/macOS: **Terminal**
2. Navigate to the Cell-ACDC folder with the command `cd` like you did when you installed it. Note that the code is in the subfolder called `src` which means you need to type an additional `cd src`
3. **Activate** the environment:
  - Conda: `conda activate acdc`
  - pip on Windows: `.\env\Scripts\activate`
  - pip on Unix/macOS: `source env/bin/activate`
4. Run the main launcher:
  - Windows: `python main.py`
  - Unix/macOS: `python3 main.py`
5. If you get the error `ImportError: No module named 'Tkinter'` you need to install tkinter with the command `sudo apt-get install python3-tk`
6. The first time, it will take 1 or 2 minutes to launch. The next times it will be faster. Once launched, you should get the following Welcome Guide window.

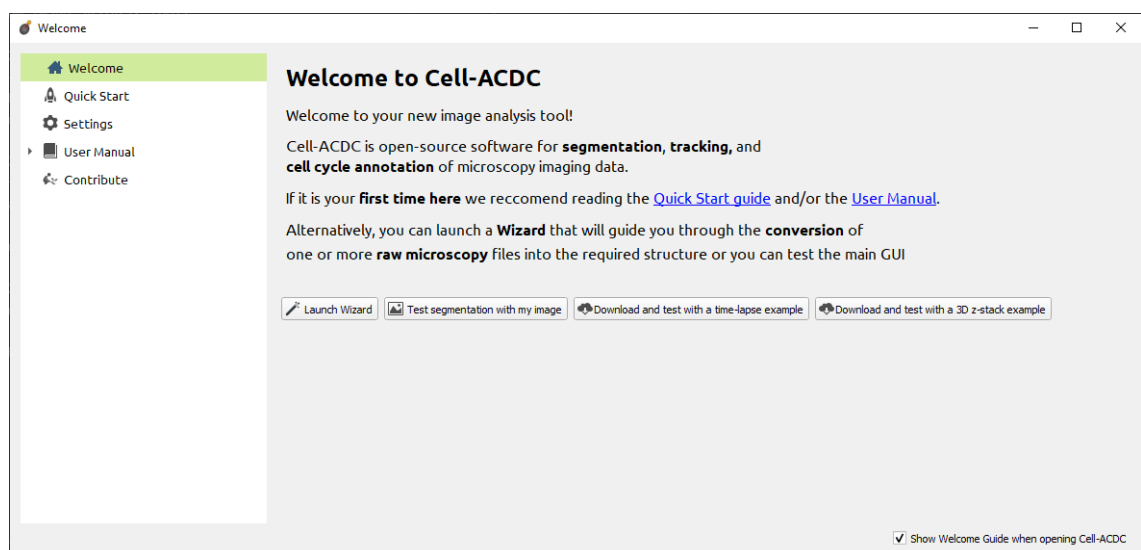




## Create required data structure from microscopy file(s)

**NOTE:** Unfortunately, loading a microscopy file to create the required data structure using Cell-ACDC is functional **only on Windows OS**. We are working on extending support to other Operating Systems. If you are not on Windows check out the next section [“Create data structure using Fiji Macros”](#)

To load a microscopy file, Cell-ACDC uses the java library [Bio-Formats](#) and the Python library [python-bioformats](#). The python-bioformats library was developed for [CellProfiler](#) and it is **embedded** into Cell-ACDC. It essentially allows you to run the java code from Bio-Formats from Python. Have a look [here](#) for a list of supported file formats.

To load a microscopy file into the Cell-ACDC pipeline, we first have to **convert it into a specific data structure**. We included a module that allows you to **automatically create the required data structure**. However, if it fails, you can create it manually with ImageJ/Fiji. Read the section of this manual called [“Manually create data structure from microscopy file\(s\)”](#).

1. From the [main launcher](#) (could be behind Welcome Guide window) click on the “Create data structure from microscopy file(s)...” button and follow the instructions of the Wizard.

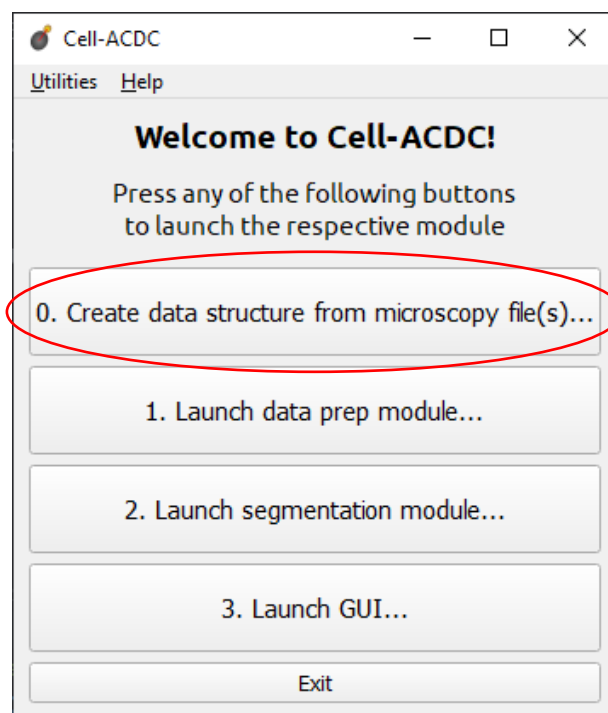

2. Once the creation of the data structure is finished, you are ready to start using your new labelling tool! The easiest way to start is from the [Quick Start](#) section of this User Manual.

## Create data structure using Fiji Macros


1. If you don't have Fiji, download and install it from [here](#)
2. Open the Fiji app. Go to File → Open... and then navigate to the `/Cell_ACDC/FijiMacros` folder.
  - a. If you have a **single microscopy file for each experiment** containing multiple positions (with or without multiple channels) open the `single_file.ijm` file.
  - b. If you have **multiple microscopy files for each experiment**, (one file for each position/series) open the `multiple_files.ijm` file. In this case, you need to put all the microscopy files from the same experiment into an empty folder.
3. Before running the macro, you have to **edit the channels**. At about line 3 modify the `channels` variable by writing a name for each channel in your file. They **MUST be in the exact same order they are in the original microscopy file**. If you don't know the order, open the file in Fiji first and check the order of the channels there.
4. Run the script with the `Run` button. The script will ask you to select a folder containing the microscopy files (`multiple_files.ijm`) or a microscopy file (`single_file.ijm`). Select the folder/file and wait. You can see the progress in the Log Window. The script will silently open the microscopy files and save each Position (if the file contains multiple positions) into a separate folder.

If successful, you should now have a **new folder for each position** called Position\_1, Position\_2 etc. Inside each Position folder you should have a folder called **Images**. Inside Images you now have **one .tif file for each channel**.


Check that each channel name (appended at the end of the .tif file) corresponds to the actual data contained in the .tif file. If it doesn't correspond, delete all the created folders, and repeat the process with the **right order of channel names** in the variable called `channels` (point 3.).

## Manually create data structure from microscopy file(s)

6. If you don't have Fiji, download and install it from [here](#)
7. Open the Fiji app and launch the Bio-Formats importer from the menu "Plugins → Bio-Formats → Bio-Formats Importer" and select your microscopy file (one at the time).
8. Check the options "Use virtual stack" and "Split channels" as in the screenshot below

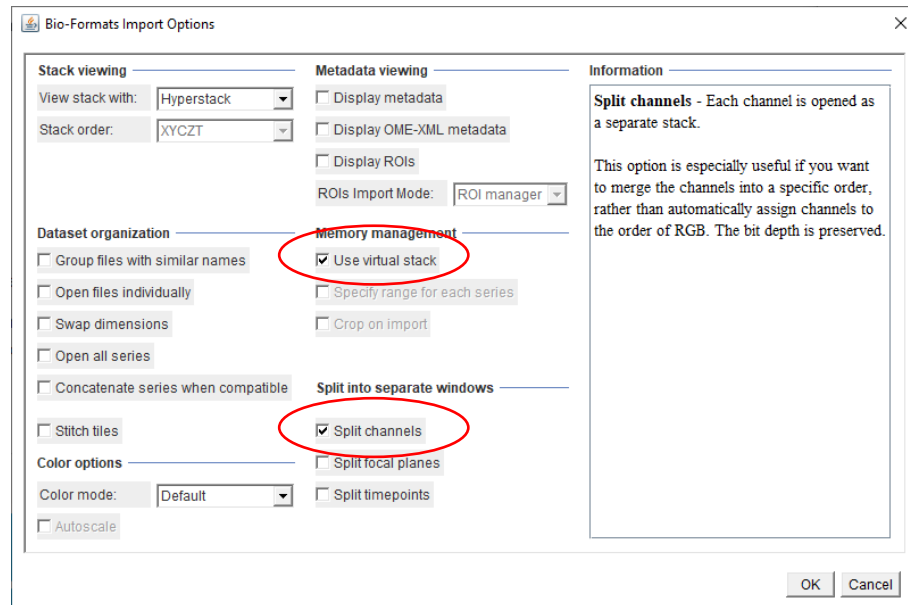

9. If you have multiple positions (series) in the file you opened, you will be asked to select a position. We recommend opening one position at the time, to avoid memory issues.
10. You should now have one image window for each channel you had in the file. Select the window with the image data from the first channel (the window name should be something like "filename ... C=0"), then "File → Save as... → Tiff"
11. As a filename we recommend calling it with the same name of the original microscopy file (if it is not too long) **WITHOUT the extension** plus something like "\_s01\_channel0.tif", where instead of "channel0" you can write whatever you like (e.g., DAPI or GFP etc.) and "s01" is for first position. So, for example the phase contrast channel of a .czi (Zeiss microscope) file called ASY15-1\_15nM-01.czi can be save as ASY15-1\_15nM-01\_s01\_phase\_contr.tif
12. Save the .tif file to a path called "/Position\_1/Images".
13. Repeat 2-7 for all the other positions.
14. In the end you should have the following folder structure:

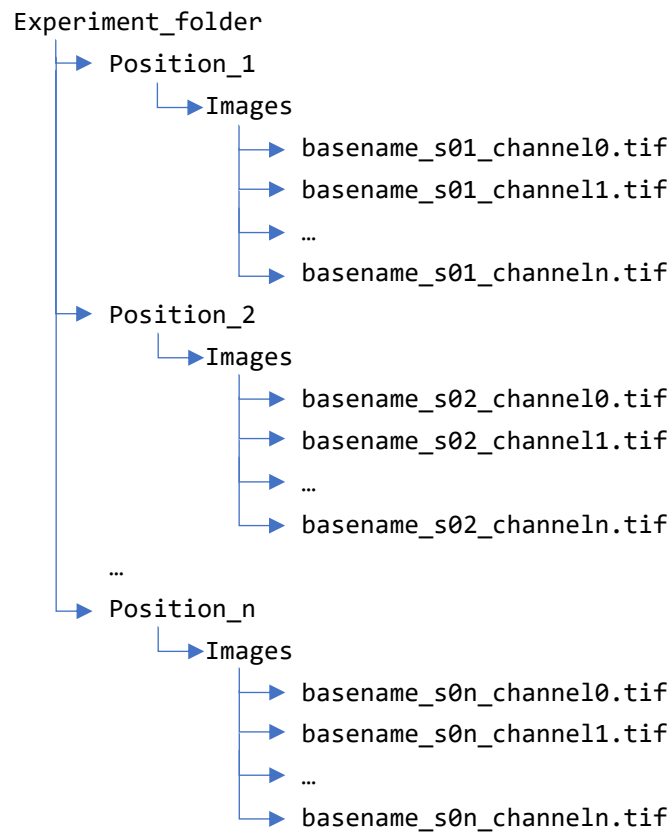

15. Once the creation of the data structure is finished, you are ready to start using your new labelling tool! The easiest way to start is from the [Quick Start](#) section of this User Manual.


## Quick Start

Cell-ACDC is composed of **three main modules**:

- **Data prep**: align time-lapse data, crop, and select a z-slice or z-projection and a ROI for segmentation. [More details...](#)
- **Segmentation**: to automatically segment multiple experiments and multiple positions with the embedded deep learning models ([YeaZ](#) for yeast cells and [Cellpose](#) for various model organisms). [More details...](#)
- **Main GUI**: to visualize segmentation masks, correct segmentation and tracking errors, and cell cycle annotations. [More details...](#)

The easiest way to start is to **open the main GUI**.

Next, if you already created the data structure (see [Load microscopy file](#) section) you can click on the “Open Folder” button on the toolbar, otherwise you can go to “File → Open image/video file...”. To start the main GUI, click the “Launch GUI...” button on the main launcher.

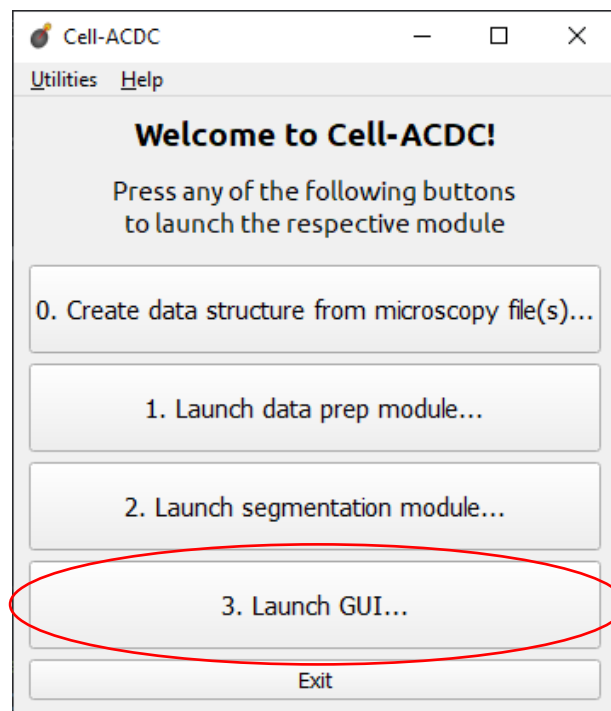

# Data Prep module


To use the data prep module, you need to **first create the required data structure**. See [this](#) section.

Use the **Data Prep** module if you need to do one of the following **tasks**:

- a) Select a **z-slice** or **z-projection** for segmentation of 3D z-stacks.
- b) **Align frames** of time-lapse microscopy data (RECOMMENDED, it is revertible).
- c) Calculate **background metrics** (median, mean etc.) from one or more **rectangular areas**. The median will be used later for background subtraction. The areas are movable and resizable.
- d) Select a region of interest (**ROI**) for segmentation.
- e) **Crop** images to reduce memory usage (RECOMMENDED, if possible).

## Loading data

1. **Launch** the data prep module, click on the “1. Launch data prep module...” button on the main launcher.
2. Click on the “**Open Folder**” button on the toolbar.

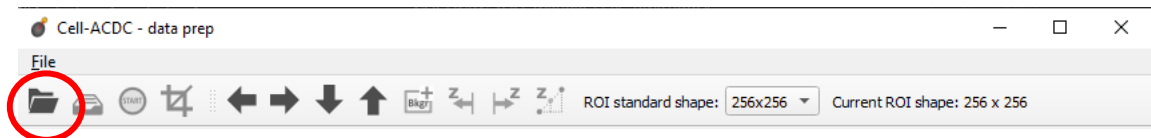

3. Select a specific Position folder or the entire experiment folder.
4. Follow the instructions in the pop-up windows. Make sure to enter the **correct metadata**.

*NOTE: For time-lapse microscopy you can load only one position at a time. Select multiple positions only if you have single 3D z-stacks or single 2D images.*


## Usage

1. If your data does not contain 3D z-stacks go to point 2. Otherwise, you can visualize the **z-slices** with the **scrollbar** below the image or choose a **z-projection** method with the selector on the right side of the scrollbar.


Every time you change the visualization method, the system will save it. It will then assume that the **last visualization is the preferred one** and it will be used for **segmentation**.


For **time-lapse data** you have **additional buttons** to help with the selection. Go to the section “Additional functions” for details about their functionality.

- If you do not need to select a ROI, crop, align or calculate background metrics you can close the window. Otherwise press the **“Start”** button on the toolbar and follow the instructions on the pop-up windows.

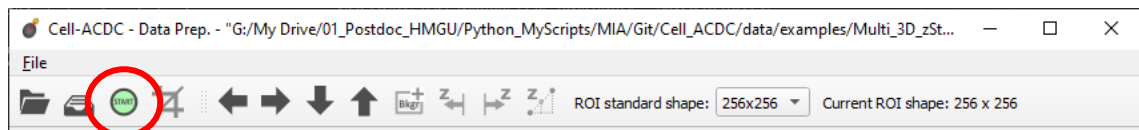

- The GUI now will be **unresponsive** until the process terminates, so do not close it. You can follow **progress** in the **terminal**. Once it finishes, a **red rectangle** will appear, along with a grey rectangle (see screenshot below). If you do not need to select a **ROI**, calculate **background metrics**, or crop you can close the window now, otherwise go to the next point.

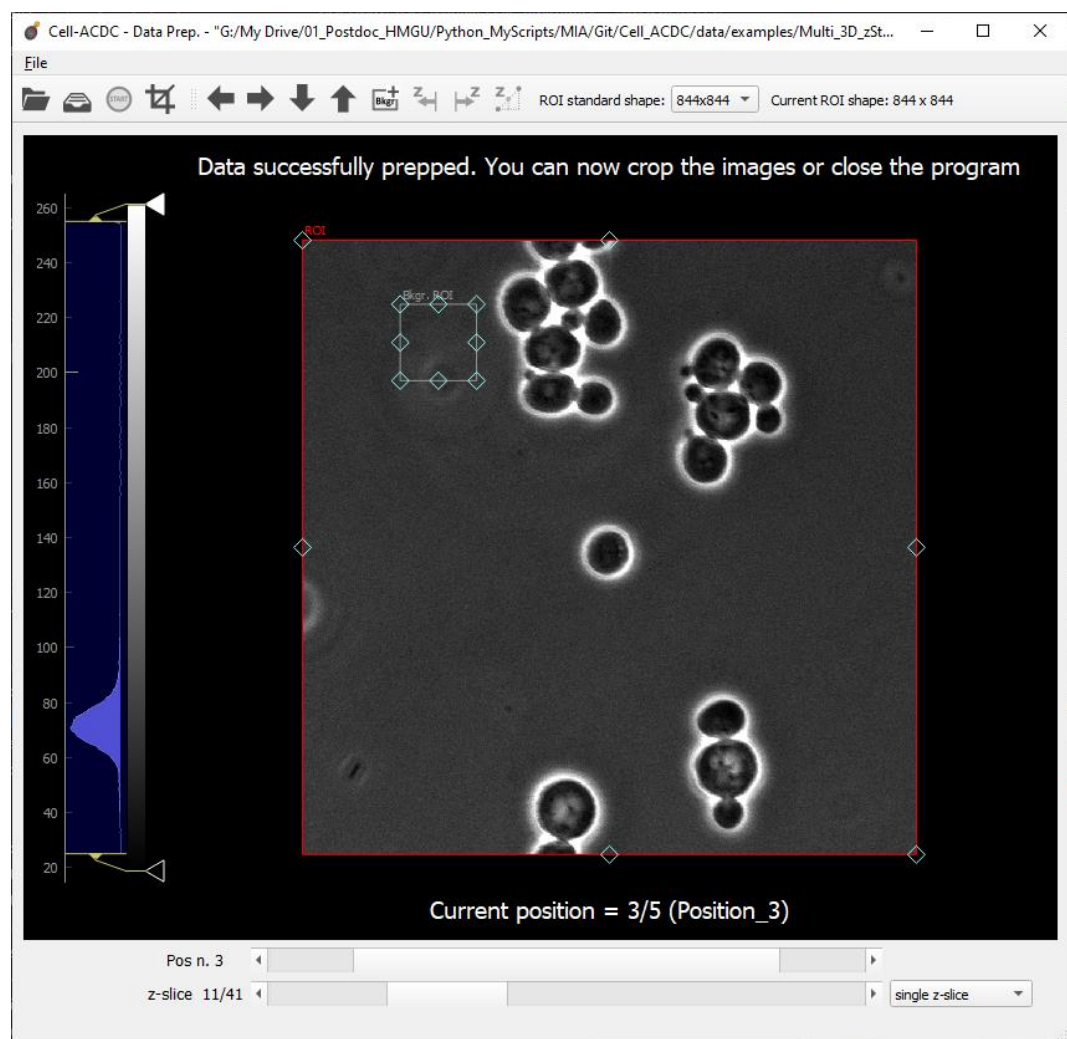

- The **red ROI** is used for either **cropping** or saving the coordinates where to compute **segmentation**. The **grey ROI** (Bkgr. ROI) is used to calculate **background metrics** from that area (median, mean, quantiles etc.). You can add more background ROIs with the “Add ROI where to calculate background intensity” () button on the toolbar.
- Resize** and **move** the ROIs until you are happy with their position and size, click on the **green tick** button on the toolbar, then follow the instructions in the pop-up windows. The GUI will be **unresponsive** until the process terminates, so **do not close it**. You can check progress in the terminal.

## Additional functions



Go to **previous/next** position or frame (time-point).



Go **10** positions or frames **backward/forward**.



Use the same **z-slice** from current frame to all past/future frames.



Use **linearly interpolated z-slices** from first frame to current frame.

ROI standard shape:

Select one of the **standard shapes** for the red ROI.



## Segmentation module


The segmentation module is used for **automatically segmenting multiple experiments and multiple positions** in one session.

To use the segmentation module, you need to **first create the required data structure**. See [this](#) section.

*NOTE: if you are just testing, you can also segment in the main GUI. Use this module when you need to segment many experiments and/or many positions.*

### Usage

To use this module, simply follow the **instructions** in the **pop-up windows**. To **launch** the module, click on the “2. Launch segmentation module...” button on the main launcher.

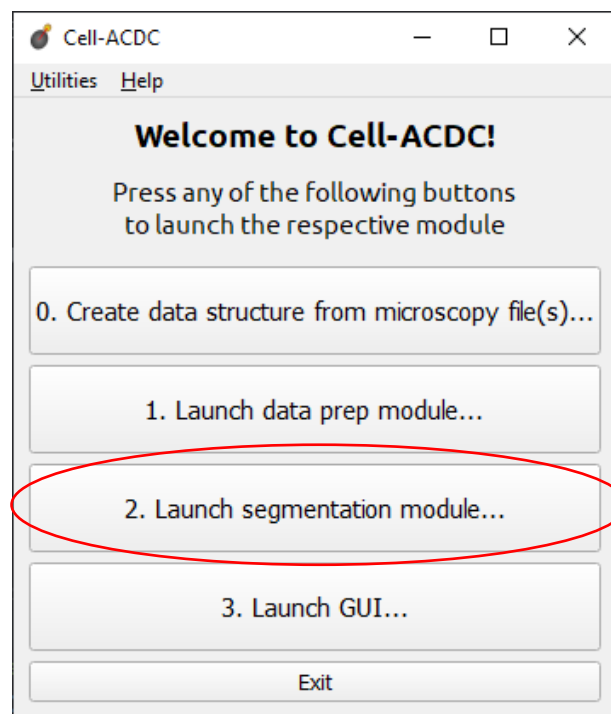

## Main GUI

The main GUI is the actual **core of Cell-ACDC**. It serves multiple purposes:

- a) **Test** which **segmentation method** works best for your dataset.
- b) **Correct** segmentation and tracking errors.
- c) Cell cycle **annotations**.

As for all the other modules, you can load **2D**, **3D** (either single z-stacks or 2D images over time) and **4D** (3D z-stacks over time) images or videos.


### Loading data

1. **Launch** the GUI module by clicking on the **3. Launch GUI...** button on the main launcher.
2. Depending on the data structure, do one of the following actions:

- If you already created the **data structure** following the instructions in [this](#) section (recommended) then click on the Open Folder button.



- If you have a single image (.tif, .png, .jpg, etc.) or video (.mov, .avi) go to **File → Open image/video file...**

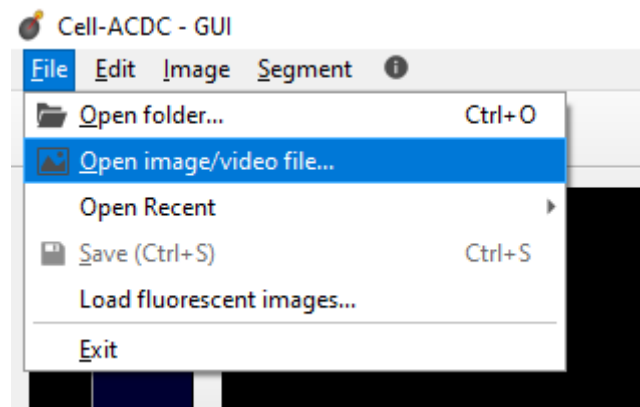

3. Select a specific file, a Position folder, or the entire experiment folder.
4. Follow the instructions in the pop-up windows. Make sure to enter the **correct metadata**.

*NOTE: if you load a **single image** or **video** file **without the required data structure**, the Cell-ACDC output will be saved in a sub-folder called **<timestamp>\_acdc\_output***

## Usage with time-lapse data

For **time-lapse data**, you can load one position (one video) at a time. With this data, the GUI has **three modes** that can be toggled from the selector on the toolbar:

- a) **Viewer** mode (default mode, used only for visualisation).
- b) **Cell cycle analysis** mode.
- c) **Segmentation and tracking** mode.

The **main idea** is that when you visit a frame for the first time, some automatic functions are triggered: **tracking** in **Segmentation and tracking** mode, **mother-bud pairing** in **Cell cycle analysis** mode.

These functions are **not triggered** when you visualize a frame that you **already visited before**. You can always call any function manually (see [this](#), [this](#), or [this](#) section of the manual).


Give a **quick read** to

Tips and tricks and then **start using the GUI**.

If you are unsure about the function of any item in the GUI, see [this](#), [this](#), or [this](#) section of the manual.

## Usage with snapshot data (no time-lapse)

For **snapshot data**, you can load **multiple positions** at the same time. When prompted, simply click on multiple selection button, and then select the positions with `Ctrl+click` for selecting specific positions, or `Shift+click` to select a range, or `Ctrl+A` to select all.

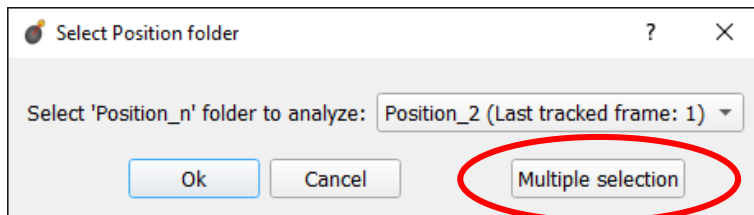

Once loaded, you can **navigate** through positions with **left and right arrow** or with the Position **scrollbar** below the left image.

Give a **quick read** to [Tips and Tricks](#) and then **start using the GUI**.


If you are unsure about the function of any item in the GUI, see [this](#), [this](#), or [this](#) section of the manual.



## Tips and tricks

- Most of the **functions** are available from the **toolbar** on the top of the window:



- Activate the function with a **SINGLE-click** on the button
- When you **hover a button** with the mouse cursor you get a **tool tip** on how to use that function:

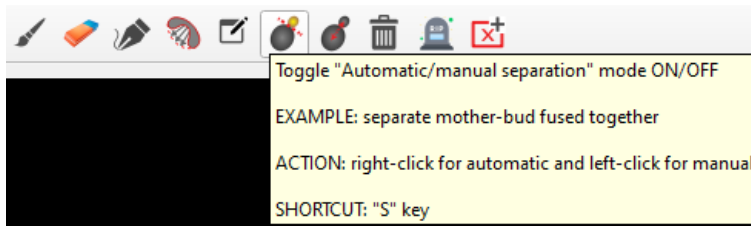

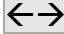- The **tool tip** will tell you whether you need **RIGHT-click** or **LEFT-click** for that function
- Functions **NOT** present on the **toolbar**:
  - *Middle-click* (scrolling wheel) → delete the segmented object you click on
  - **H** key → automatic zoom on the segmented objects
  - **Ctrl+P** → visualize cell cycle annotations in a table
  - **L** key → relabel object IDs sequentially (1,2,3...etc)
- To **test** the available **segmentation models**, use the **Segment** menu.
- To **navigate frames** (time-lapse data) or **positions** (snapshots data), use the  **arrows** on the keyboard.
- To visualize the frames of time-lapse data in a second window click on the Slideshow button on the toolbar: 
- **Personalize** settings such as **Font Size**, **overlay colour** and **text's colour** from the **Edit** menu.

## Functions activated from the toolbar

### File toolbar



**Open folder:** used to **automatically load** single or multiple positions from the standard data structure. See [this](#) section for details on how to create the required data structure.



**Save:** used to save **all data** (segmentation mask, cell cycle annotations and metrics calculated from the loaded fluorescent images).



**Show in Explorer/Finder:** open the loaded folder into your Explorer/Finder.



**Reload segmentation file:** used to reload the segmentation labels from the hard drive. Use it if you want to go back to the saved state.



**Undo (Ctrl+Z):** almost all the performed actions are undoable. Currently you can undo up to 5 actions in the past.



**Redo (Ctrl+Y):** repeat an undone action.

### Visualize toolbar




**Go to previous** frame (time-lapse data) or Position (snapshot data).  
*SHORTCUT: Left arrow on the keyboard*



**Go to next** frame (time-lapse data) or Position (snapshot data).  
*SHORTCUT: Right arrow on the keyboard*



**Open** the images in a **second window**. This window will have no annotations, which means it is the fastest way possible to visualize the frames. Handy in many situations.

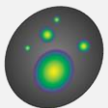

**Overlay a second signal:** when you press this button for the first time, you will be asked to choose which signal to overlay. The next times it is used to toggle **overlay on/off**.



**Ruler:** draw a line with the left button to **measure the distance** between two points (*pixels* and  $\mu m$ ). Distance will be displayed in the bottom-left corner of the window.

## Edit toolbar

**Brush** (*left-click motion*): used to modify a segmented object or paint a new object with a circular brush. Increase/decrease the size of the brush with up/down arrows on the keyboard. *SHORTCUT: B key*




Two **modes**:

- To **draw a new object** start painting on the background (brush cursor is white).
- **Modify an object** by starting to paint from the object (brush cursor takes the color of the object you are about to modify)

*NOTE: The brush will add the new object UNDER existing objects unless you first press the B key twice. The brush button will then turn red and new objects will be added ABOVE existing objects. Restore default behaviour by pressing B twice again.*

**Eraser** (*right-click motion*): used to erase parts of a segmented object with a circular cursor. *SHORTCUT: X key*

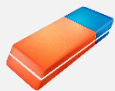

Two **modes**:

- Start erasing from the background to erase **all the objects** you will pass over (eraser cursor is red).
- Start erasing from a specific object to erase **ONLY** that object (eraser cursor keeps the colour of the object you start erasing from).

*NOTE: To enforce erasing any object you pass over, even if you start from a specific one, press the X key twice. The eraser button will then turn RED. Restore default behaviour by pressing X twice again.*

**Curvature tool**: used to draw new objects by drawing a spline with multiple anchor points. *SHORTCUT: C key*



Two **modes**:

- Consecutive **left-clicks** to **manually** draw an object with a spline.
- **Right-click** with a drawing motion to **automatically** follow an intensity line.



**Hull contour**: replace object with its hull contour image. Useful to **fill holes** and **cracks**. *SHORTCUT: F key*

**Usage**: first *activate the button* (left-click or shortcut) and then *RIGHT-click* on the requested object.



**Edit ID**: replace the ID of an object with a manually inserted one. *SHORTCUT: N key*

**Usage**: first *activate the button* (left-click or shortcut) and then *RIGHT-click* on the requested object.

**Separate objects:** used to separate merged objects. *SHORTCUT: S key*



**Usage:** first *activate the button* (left-click or shortcut) and then *RIGHT-click* on the requested object to first attempt **automatic separation** (works well with two objects separated by a constriction) or *LEFT-click* to go straight to **manual separation**.

When **manual separation** is triggered, a window with only the object you clicked one will appear (see screenshot below). To separate the object along a curved line, you need **three RIGHT clicks**:

1. First *RIGHT-click* where to **start** a straight **line**.
2. Second *RIGHT-click* where to **end** the straight **line**.
3. Third *RIGHT-click* to **set the curvature** of the separating curve.

The object will then be **separated** along the drawn line.

*NOTE: to help with deciding where to draw the line you can **overlay** the intensity signal with the overlay button*

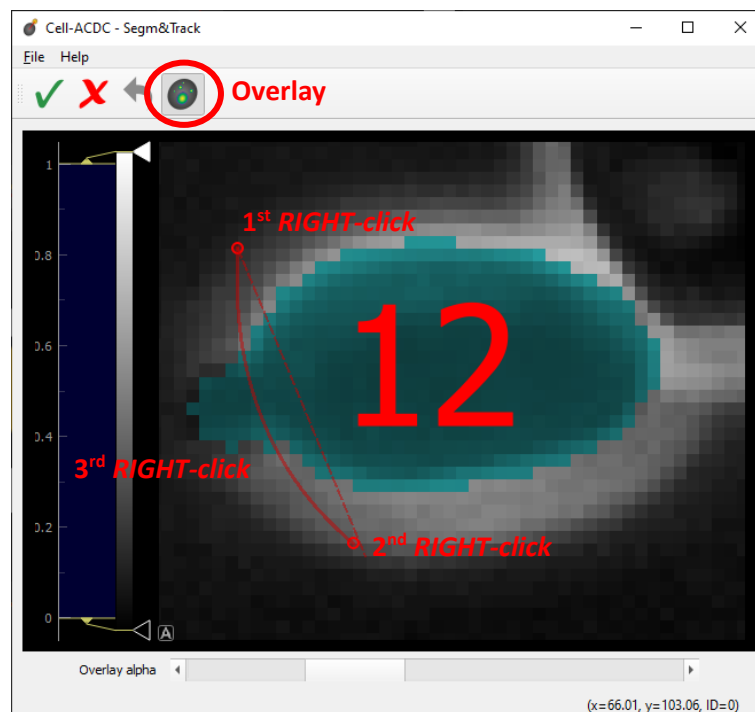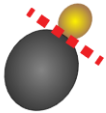

**Merge objects.** *SHORTCUT: M key*

**Usage:** first *activate the button* (left-click or shortcut) and then *drag-and-drop* motion with *RIGHT-button* between the two objects.



**Annotate cell as “Excluded from the analysis”:** the column called `is_cell_excluded` (on the saved data, see [this](#) section) will have a TRUE for the clicked object from current frame to the end. *SHORTCUT: R key*

**Usage:** first *activate the button* (left-click or shortcut) and then *RIGHT-click* on the requested object. Repeat to undo from that frame onwards.





**Annotate cell as “Dead”:** the column called `is_cell_dead` (on the saved data, see [this](#) section) will have a TRUE for the clicked object from current frame to the end. **SHORTCUT:** *R* key

**Usage:** first *activate the button* (left-click or shortcut) and then *RIGHT-click* on the requested object. Repeat to undo from that frame onwards.

**Add a delROI:** used to add a rectangular area to automatically delete all objects contained in or touched by it.



**Usage:** click on the button to add a ROI. **Resize and move** the ROI to either **delete** objects or **restore** objects when they are not contained/touched anymore. Delete the ROI with right-click on it → Remove ROI.

**NOTE:** Objects are **permanently deleted** only when you save AND close the GUI.



**Repeat tracking:** used to repeat tracking on the current frame.

**Note:** The tracking in the GUI is **different** from the tracking of the segmentation module. In the segmentation module we kept the algorithm proposed by the YeaZ model.

## Cell cycle annotations toolbar



**Assign bud/sister to mother/sister.** **SHORTCUT:** *A* key

**Usage:** first *activate the button* (left-click or shortcut) and then *drag-and-drop* motion with *RIGHT-button*. Right-click on bud and release on mother, or **right-click on bud** then **right-click on mother**.



**Annotate** that a cell does **not** have a **fully known history**. **SHORTCUT:** *U* key

**Usage:** first *activate the button* (left-click or shortcut) and then *RIGHT-click* on the requested object.

**Tip:** two examples of cells with unknown history are cells already present at the first frame and cells appearing from outside the field of view



**Reinitialize cell cycle annotations** to default. Default is all cells in **G1** without a relative assigned to it.

## Functions activated from the menus

File → Load fluorescent images...

Used to **load** as many **additional images** (e.g., fluorescence signal) as you want.

Loaded images will be used to calculate metrics such as mean, median, total amount etc. See [this](#) section for more details.

Edit → Smart handling of enabling/disabling tracking

The GUI has built-in **automatic tracking** for **time-lapse data** with the following behaviour:

- If tracking is active (Disable tracking checkbox on the toolbar is UNCHECKED) when you visit a frame that you have **never visited before**, objects will be **automatically tracked** compared to previous frame.
- When you visit a frame **already visited** before, it will **not be tracked**.

You can disable this automatic behaviour by unchecking Smart handling of enabling/disabling tracking

When you **disable** the **smart** handling, you can **enforce tracking** on all visited frames no matter if they were previously visited or not. To enforce, use the Disable tracking checkbox on the toolbar.

*Tip: useful when you know you have to repeat tracking on already visited frames.*

Image → Normalise intensities → ...

You can choose to **normalise the intensities** of the displayed images (saved data will not be modified) with the following methods:

- Do not normalise. Display **raw image**.
- Convert to **floating point** format with values [0, 1] → simply convert to floating point, no normalisation involved.
- **Rescale** to [0,1] → intensities are first converted to floating point if needed and then STRETCHED to convert the entire [0,1] range.
- **Normalize by max** value: divide by the max of the intensities.

## Additional functions

- *Middle-click* (scrolling wheel) → delete the segmented object you click on.
- H key → automatic zoom on the segmented objects.
- Ctrl+P → visualize cell cycle annotations in a table.
- L key → relabel object IDs sequentially (1,2,3...etc).

## Cell-ACDC output data

**Files** saved by Cell-ACDC for a **fully analysed experiment** (example with original raw microscopy file called `Example1` and first position `_s01_`)

*NOTE: not all files are always present, it depends on whether you have 3D data, time-lapse data, or you aligned or cropped.*

<code>Example1_s01_acdc_output.csv</code>	Main table containing <b>cell cycle annotations</b> and additional <b>metrics</b> such as mean, median etc. for all the loaded channels plus all the <b>region properties</b> (as calculated by <a href="#">skimage.measure.regionprops</a> ) for each segmented object.
<code>Example1_s01_act1.tif</code>	.tif file for the channel called <code>act1</code>
<code>Example1_s01_cdc10.tif</code>	.tif file for the channel called <code>cdc10</code>
<code>Example1_s01_last_tracked_i.txt</code>	<b>Last visited frame</b> in “Segmentation and Tracking mode” with the main GUI
<code>Example1_s01_metadata.csv</code>	Table containing the <b>metadata</b> such as number of frames, number of z-slices, pixel size etc.
<code>Example1_s01_phase_contr.tif</code>	.tif file for the channel called <code>phase_contr</code>
<code>Example1_s01_segm.npz</code>	<b>Segmentation labels.</b> This is a numpy array (compressed).
<code>Example1_s01_segmInfo.csv</code>	Table containing information such as which <b>z-slice</b> or <b>z-projection</b> was used for segmentation or saving metrics of each channel.
<code>Example1_s01_align_shift.npy</code>	Numpy array containing the <b>shifts</b> applied to each frame when <b>aligning</b> . Useful for reverting to non-aligned state.
<code>Example1_s01_dataPrepROIs_coords.csv</code>	Table containing the <b>coordinates of the ROI</b> that was used to either <b>crop</b> , or <b>segment</b> only in the ROI. This is created in the data prep or segmentation stage.
<code>Example1_s01_phase_contr_aligned.npz</code>	Aligned data for the channel called <code>phase_contr</code> . This is a numpy array (compressed).
<code>Example1_s01_phase_contr_aligned_bkgrRoiData.npz</code>	<b>Data</b> from the <b>background ROIs</b> generated at the data prep stage.
<code>Example1_s01_phase_contr_dataPrep_bkgrROIs.json</code>	<b>Coordinates</b> of the <b>background ROIs</b> generated at the data prep stage

## Adding segmentation models to the pipeline

Adding segmentation models requires few minutes:

1. Create a new folder with the model's name (e.g., YeastMate) inside the `Cell-ACDC/src/models` folder
2. Create a `__init__.py` file inside the model folder, and paste the following code in it:

```
import os
import sys

model_path = os.path.dirname(os.path.abspath(__file__))

sys.path.insert(0, model_path)
```

Add any line of code needed to initialize correct import of the model.

3. Create a new file called `acdcSegment.py` in the model folder with the following template code:

```
import module1
import module2

class Model:
    def __init__(self, **init_kwargs):
        script_path = os.path.dirname(os.path.realpath(__file__))
        weights_path = os.path.join(script_path, 'model', 'weights')

        self.model = MyModel(
            weights_path, **init_kwargs
        )

    def segment(self, image, **segment_kwargs):
        lab = self.model.eval(image, **segment_kwargs)
        return lab
```

Have a look at the already implemented models (YeaZ, Cellpose, and YeastMate) as an example.

4. Create a folder inside the model folder with the weights and/or configuration files of your model. Typically, I call this folder `model`.

That's it. Next time you launch the segmentation module (or from the GUI menu `Segment`) you will be able to select your new model.

**The model parameters will be automatically inferred** from the class you created in the `acdcSegment.py` file, and a widget with those parameters will pop-up. In this widget you can set the model parameters (or press Ok without changing anything if you want to go with default parameters).

# Downstream analysis

## Overview

We provide a notebook for downstream analysis called *cell\_cycle\_analysis.ipynb* which is meant as a starting point to generate all kinds of visualizations of the data generated by Cell-ACDC.

To **start** the notebook, open a terminal and navigate to the Cell-ACDC folder. Then type `jupyter-lab`.

This should open the project folder in the jupyter notebook interface JupyterLab. In the interface, navigate into the folder `notebooks/` and open *cell\_cycle\_analysis.ipynb*.

If we now switch to the “Table of contents” panel in the left, we get an overview over the notebook’s contents:

The screenshot shows the JupyterLab interface with the notebook *cell\_cycle\_analysis.ipynb* open. The left sidebar contains the 'Table of contents' panel, which lists the following sections:

- 0 configurations
- 1 load data and perform all needed calculations on image data
- 2 Plot gallery - timelapse data
- (Volume) growth in G1 vs. mother+daughter growth in S (1st generation)
- Volume at birth vs. G1 duration (1st generation)
- Volume at birth vs. Signal concentration at birth (1st generation)
- G1 vs. S duration (1st generation)
- Plot gallery - z-stack data, no timelapse
- 2 Distribution of Cell volumes
- Volume vs. Signal amount
- Volume vs. Signal concentration
- ACDC paper figures
- Fluorescence Signal over time (centered on bud emergence)
- Volume at birth and division vs. mCitrine amount at birth (single cell) and division (combined)

The main area of the notebook shows the following code cells:

```
[1]: import os
import sys
sys.path.append('../src/')
import glob
import numpy as np
import pandas as pd
pd.set_option("display.max_columns", 150)
pd.set_option("display.max_rows", 30)
import matplotlib.pyplot as plt
import matplotlib.patches as mpatches
import matplotlib.lines as mlines
import seaborn as sns
sns.set_theme()
import cca_functions
```

```
[2]: %load_ext autoreload
%autoreload 2
```

### configurations

- follow the file selection dialog:
  - select microscopy folder in first step
  - select positions of the selected folder in second step
- repeat to add more positions to the analysis
- positions selected within one iteration of the dialog will be pooled together in the following analyses

```
[3]: data_dirs, positions = cca_functions.configuration_dialog()
file_names = [os.path.split(path)[-1] for path in data_dirs]
image_folders = [[os.path.join(data_dir, pos_str, 'Images') for pos_str in pos_list] for pos_list, data_dir in zip(pos_list, data_dirs)]
# determine available channels based on first(!) position.
# Warn user if one or more of the channels are not available for some positions
first_pos_dir = os.path.join(data_dirs[0], positions[0][0], 'Images')
first_pos_files = os.listdir(first_pos_dir)
channels, warn = cca_functions.find_available_channels(first_pos_files, first_pos_dir)
```

```
[4]: overall_df, is_timelapse_data, is_zstack_data = cca_functions.calculate_downstream_data(
    file_names,
    image_folders,
    positions,
    channels,
    force_recalculation=False
)
```

Load files for Multi\_3D\_zStacks\_labeled, Position\_1...

Number of cells in position: 17

Number of annotated frames in position: 1

Calculate regionprops on each frame based on Segmentation...

Calculate mean signal strength for every channel and cell...

Saving calculated data for next time...

Load files for Multi\_3D\_zStacks\_labeled, Position\_2...

Number of cells in position: 12

Number of annotated frames in position: 1

Calculate regionprops on each frame based on Segmentation...


Calculate mean signal strength for every channel and cell...

Saving calculated data for next time...


Depending on your input data, you can generate plots based on [timelapse data \(1\)](#) or [zStack data \(2\)](#) which does not contain data over time. In any case, you should run the **configuration and data load step (0)** before.

## Configuring the analysis

Running the configuration cell opens a widget (**possibly in the background**) where we can select the data folder containing the “*Position\_n*” folders (see [First steps](#)):

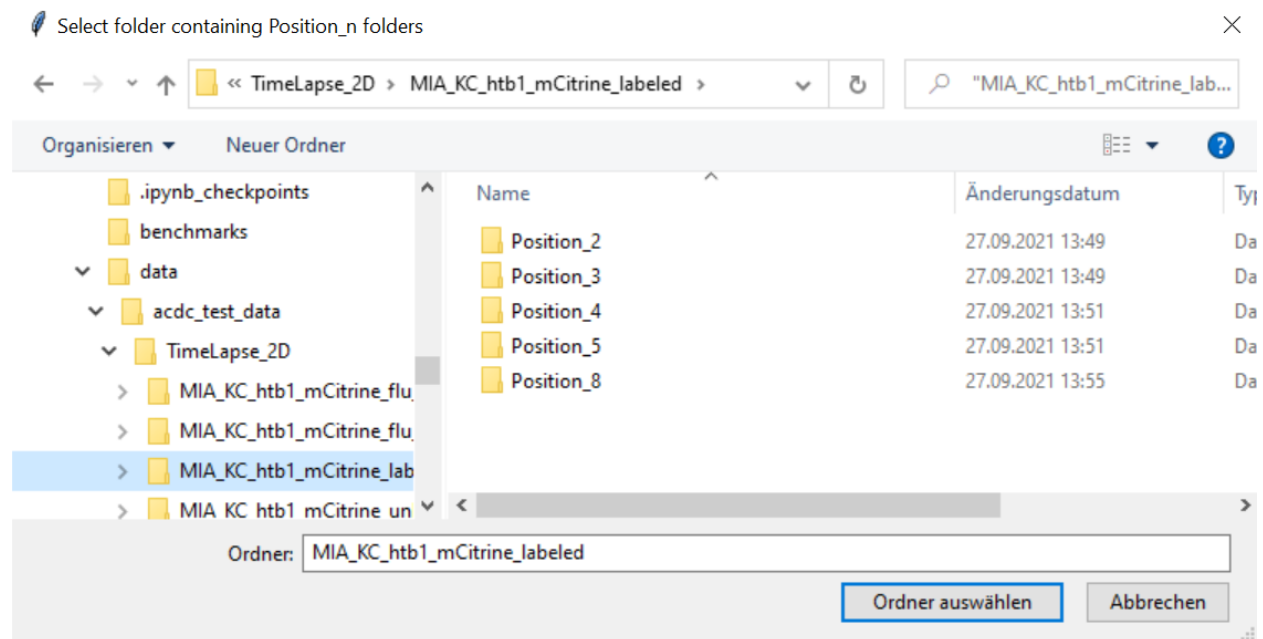


Selecting the folder leads to a window where we can select positions:

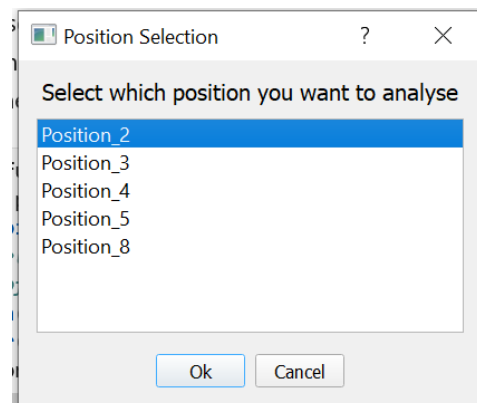

Multiple positions can be selected with **Ctrl/Shift**. All positions can be selected with **Ctrl+A**.

Next, we can select if we want to add another file and the corresponding positions to the analysis:

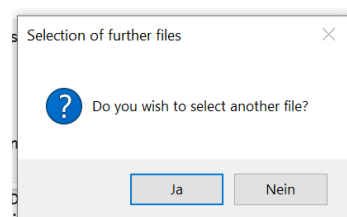

Additional files added when choosing “Yes” here, will be added to a second/third/... pool which can be helpful for comparing multiple experiments in the notebook’s plots. Choosing “No” here concludes the file selection and configuration of the notebook.

## Pre-processing

Pre-processing is performed by executing the next two cells:

```
[*]: overall_df, is_timelapse_data, is_zstack_data = cca_functions.calculate_downstream_data(
    file_names,
    image_folders,
    positions,
    channels,
    force_recalculation=True
)
```

[illegible]

```
[*]: # if cell cycle annotations were performed in ACDC, extend the dataframe by a join on each cells relative cell
      if 'cell_cycle_stage' in overall_df.columns:
          overall_df_with_rel = cca_functions.calculate_relatives_data(overall_df, channels)
      # If working with timelapse data build dataframe grouped by phases
      group_cols = [
          'Cell_ID', 'generation_num', 'cell_cycle_stage', 'relationship', 'position', 'file',
          'max_frame_pos', 'selection_subset', 'max_t'
      ]
      # calculate data grouped by phase only in the case, that timelapse data is available
      if is_timelapse_data:
          phase_grouped = cca_functions.calculate_per_phase_quantities(overall_df_with_rel, group_cols, channels)
          # append phase-grouped data to overall_df_with_rel
          overall_df_with_rel = overall_df_with_rel.merge(
              phase_grouped,
              how='left',
              on=group_cols
          )
          phase_grouped.head()
```

Data is calculated if `force_recalculation` is set to `True` or if there is no pre-calculated data to load. After calculation, the data gets saved and automatically loaded the next time the notebook is started. By default, a check is performed and data is loaded from disk if it is already available (`force_recalculation=False`).

## Plot gallery – Time-lapse data


Once pre-processing is finished, plots in the gallery can be produced by executing the corresponding cells:

### Volume at birth vs. Signal concentration at birth (1st generation)

```
[12]: # set channel name here:
ch_name = 'mCitrine'
# obtain table where one cell cycle is represented by one row:
# first set of columns (like phase_length, growth...) for G1, second set of cols for S
plot_data4 = phase_grouped[phase_grouped.cell_cycle_stage=="G1"]
plot_data4 = plot_data4[plot_data4.complete_phase==1]
plot_data4 = plot_data4[plot_data4.generation_num==1]

sns.set_theme(style="darkgrid", font_scale=2)
# Initialize the figure
g = sns.lmplot(x="phase_volume_at_beginning", y=f"phase_{ch_name}_concentration_at_beginning", data=plot_data4,
               hue="selection_subset", height=10, )
g._legend.set_title('Position Pool')
g.set(yscale="log")
ax = plt.gca()
ax.set_ylabel("mCitrine signal amount per volume in cell [a.u.]", fontsize=20)
ax.set_xlabel("Volume at birth (first cytokinesis) [fL]", fontsize=20)
ax.set_title("Volume at birth vs mCitrine signal amount per volume (1st generation)", fontsize=30)
plt.show()
```

### Volume at birth vs mCitrine signal amount per volume (1st generation)

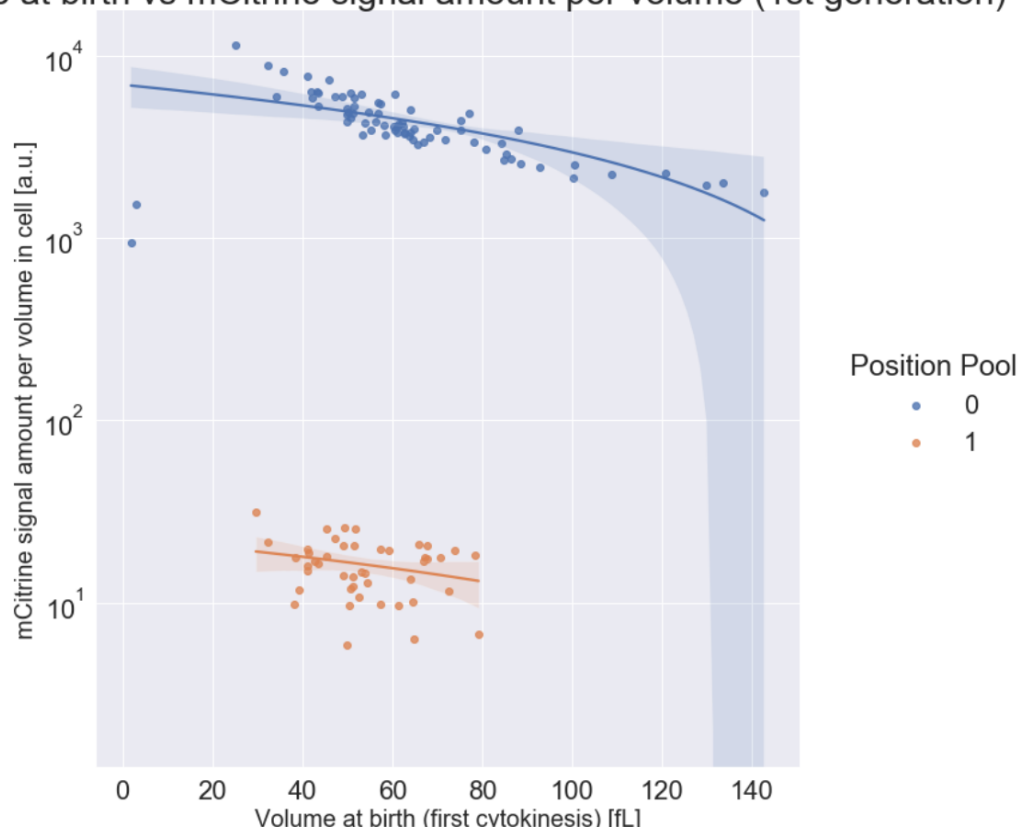

This example shows the *Htb1*-mCitrine amount per volume on the y axis and the volume at birth on the x-axis. Clearly the tagged strain (position pool 0) shows higher signal amounts as the auto-fluorescent control showing almost no signal.




## Plot gallery – 3D z-stacks data

Also 3D stacks can be loaded into the downstream analysis notebook and visualized in various ways:

```
[9]: sns.set_theme(style="ticks", font_scale=2)

# Initialize the figure
plt.figure(figsize=(10,10))
sns.histplot(
    x='cell_vol_fl',
    data=overall_df,
    hue='relationship',
    bins=20,
    legend=False
)
ax = plt.gca()
labels = [
    'Mother cells',
    'Buds'
]
handles = [
    mpatches.Patch(color=sns.color_palette('pastel')[0]),
    mpatches.Patch(color=sns.color_palette('pastel')[1])
]
ax.legend(
    handles=handles,
    labels=labels,
    loc='upper right',
    #bbox_to_anchor = (1,0.2),
    framealpha=0.5
)

# Tweak the visual presentation
ax = plt.gca()
ax.set_xlabel("Cell volume [fL]", fontsize=20)
ax.set_title(f"Volume distribution, n: {overall_df.shape[0]}", fontsize=30)
#sns.despine(trim=True, Left=True)
plt.show()
```

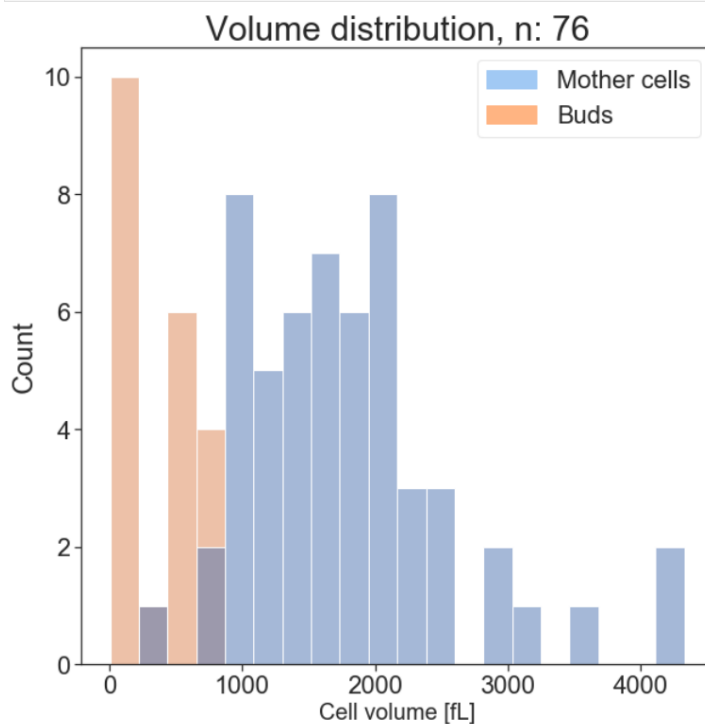

Unsurprisingly, the buds in the test position show way smaller volumes than mother cells. Outliers like the mother cell showing a very small volume in this case, can be identified and checked in the GUI.